











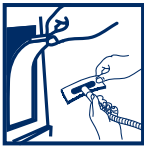






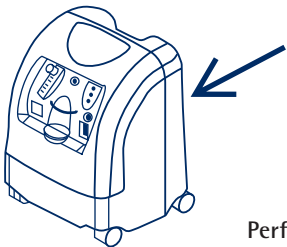
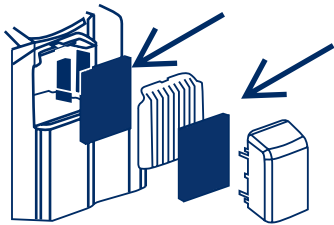
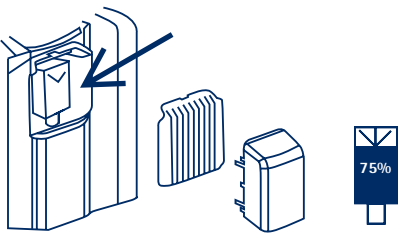
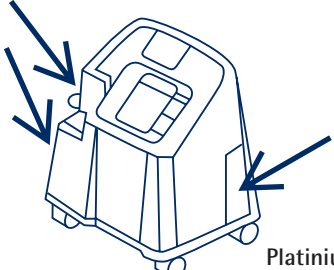
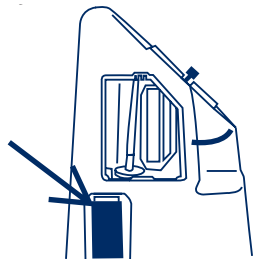
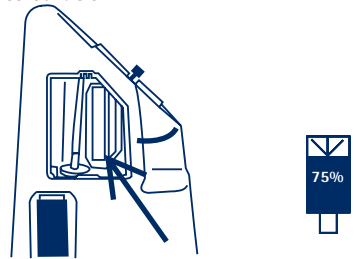
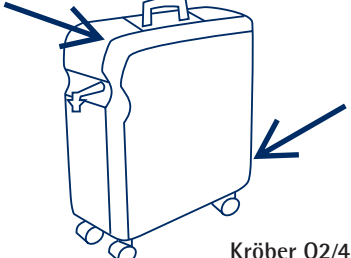

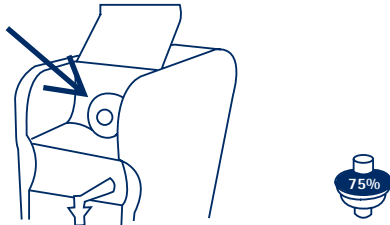
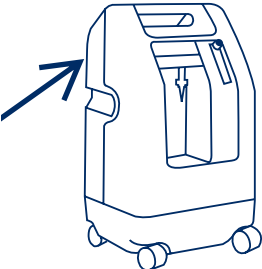
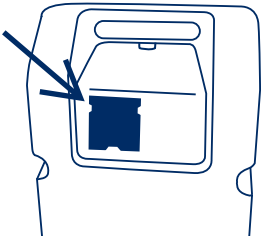
Westfalen Medical BV
Rigastraat 14
7418 EW Deventer
Tel. 0570 - 858 450
Fax 0570 - 858 451
www.westfalenmedical.com
info@westfalenmedical.nl

Instruction card for oxygen concentrator

<p>1. Always read the instructions.</p> 	<p>2. Do not smoke. No open fire.</p> 	<p>3. Ventilate the room well.</p> 
<p>4. Do not place the concentrator under the coat rack or against the curtains. This is to prevent oxygen saturation.</p> 	<p>5. Place the concentrator in a safe place as indicated. Always use the concentrator in the upright position.</p> 	<p>6. Do not use oil/grease.</p> 
<p>7. Do not use any electrical equipment.</p> 	<p>8. Maximum length connection hose.</p> 	<p>9. We recommend replacing the nozzle once a week.</p> 
<p>10. Turning on the concentrator. Stick the plug into the socket and turn on the concentrator.</p> 	<p>11. Configure the flow meter as prescribed by your physician.</p> 	<p>12. Check the flow by placing the nozzle and/or hose in a glass of water. In the absence of bubbles, inspect the hose and flow meter. In the absence of flow, contact Westfalen Medical.</p> 
<p>13. The coarse filter should be cleaned weekly. It can be found on the side and/or back of the concentrator.</p> 	<p>14. Do not use any extension cords.</p> 	<p>15. Turn off the concentrator after use. Do not run the concentrator when the device is not in use.</p> 
<p>16. Technical maintenance must be performed by our service staff only.</p> 	<p>17. Defects and/or questions:</p> <p>Tel. 0570 - 858 450</p> 	<p>18. Malfunctions and/or questions:</p>

Westfalen Medical BV
Rigastraat 14
7418 EW Deventer
Tel. 0570 - 858 450
Fax 0570 - 858 451
www.westfalenmedical.com
info@westfalenmedical.nl

Instruction chart for filters

<p>1. Filters are located at the back.</p>  <p>Perfect O2</p>	<p>2. Vacuum the coarse filter weekly.</p> 	<p>3. Replace the intake filter behind the coarse filter at 75% black discolouration.</p> 
<p>1. Filters are located on both sides.</p>  <p>Platinum S/9</p>	<p>2. Vacuum coarse filters on both sides weekly.</p> 	<p>3. Replace the intake filter at 75% black discolouration.</p> 
<p>1. Filters are located on both sides.</p>  <p>Kröber O2/4.0</p>	<p>2. Vacuum the coarse filter weekly.</p> 	<p>3. Replace the intake filter at 75% black discolouration.</p> 
<p>1. Filters are located at the back.</p>  <p>DeVilbiss</p>	<p>2. Vacuum the coarse filter weekly.</p> 	<p>3. Replace the intake filter at 75% black discolouration.</p> 